

Human Cryopreserved Hepatocytes

Lot number: HUM4192

Date (D/M/Y): 25/05/2017


Lonza

Triangle Research Labs – Now Part of Lonza

Cryo Characterization Report (CCR)

Lot Overview		
Qualification	Catalog Number	Manufacture Date (D/M/Y)
Cryopreserved human hepatocytes, Induction	HUCPI	04-08-2016

Storage Conditions: <-150°C


Donor Demographics								
Sex	Race	Age	BMI	Tobacco Use	Heavy Alcohol Use	Drug Use	Serological Data	Cause of Death
Female	Asian	16 yrs	15.6	No	No	No	Negative	Anoxia due to developmental complications

Additional donor demographic information, including relevant medical and medication history, is available upon request


Post-thaw Viability and Cell Quality Assessment			
Thawing Medium Used	Optimal Centrifuge Conditions	% Viability (post-thaw)	Viable cell yield per vial
MCHT50	100 g x 8 minutes	89	11 x 10 ⁶

Monolayer Assessment					
Plating Medium Used	Well Format	Culture Medium Used	Optimal Seeding Density	Initial Attachment Efficiency	Monolayer Confluency @ 120 hrs
Hepatocyte Plating Medium (MP100)	24 well	Hepatocyte Maintenance Medium (MM250)	0.8x10 ⁶ /mL	85%	100%

Characterization was completed in a 24-well collagen-coated plate with extra-cellular matrix overlay

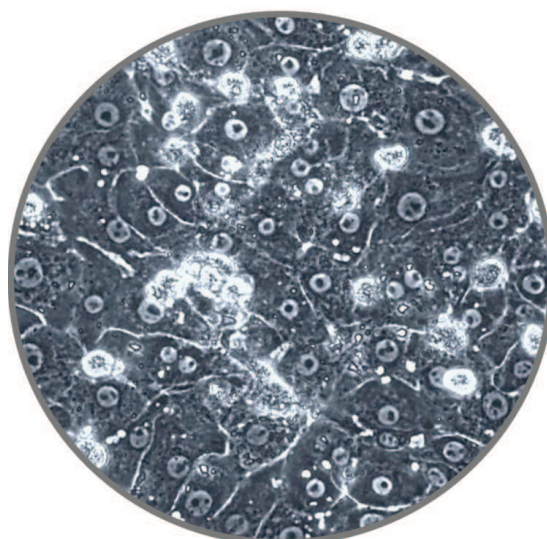

HUM4192, Post-Thaw 48 hrs, 20X

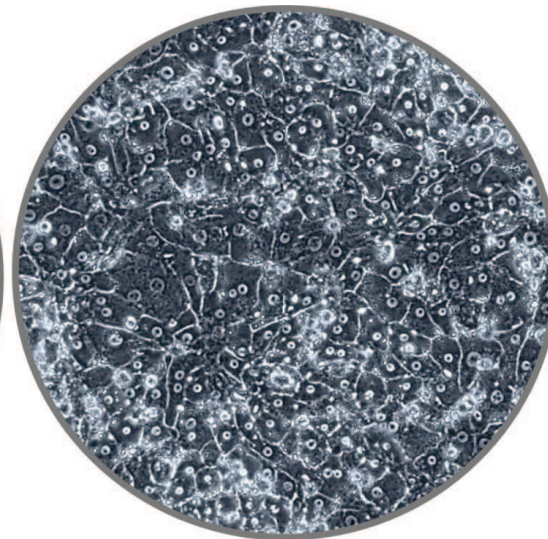

HUM4192, Post-Thaw 96 hrs, 10X

<i>Induction</i>				
Isoforms	Control Inducer	Fold Induction Specific Activity	Fold Induction mRNA Expression	Basal Induction (pmol/million cells/min)
<i>CYP1A2</i>	<i>50µM Omeprazole</i>	<i>56.1</i>	<i>77</i>	<i>21.5</i>
<i>CYP2B6</i>	<i>1mM Phenobarbital</i>	<i>32.2</i>	<i>31</i>	<i>19.6</i>
<i>CYP3A4</i>	<i>10µM Rifampicin</i>	<i>67.4</i>	<i>27.9</i>	<i>67.0</i>

Cryopreserved human hepatocytes were thawed and plated on 24-well collagen I coated plates, overlaid with Matrigel®, then dosed in triplicate with vehicle control (0.1% DMSO) or control inducers for 72 hours. The fold induction was calculated by dividing the induced level by the vehicle control level.

Media products used for characterization:

MCHT50 - Cryopreserved human hepatocyte thawing medium, 50mL

MCHT50P - Cryopreserved pooled human hepatocyte thawing medium, 50mL

Contact customer service to place an order or to obtain additional information on any of our lots. This may include supplementary donor demographic information, current inventory, and photomicrographs at multiple timepoints and magnifications.

To contact Lonza:

Customer Service: +1 800-638-8174 | order.us@lonza.com

Scientific Support: +1 800-521-0390 | scientific.support@lonza.com | Web: www.lonza.com

EU/Row Scientific Support: +32 87 321 611 | scientific.support.eu@lonza.com;

Lonza.com/TRLheps